

**PATHFINDER**





- High-quality stainless steel.
- Perfect for initial negotiation or for opening blocked canals.
- Pathfinder is sized between a #6 and #8 hand file.
- Fewer files are needed – less buckling and more efficient cutting.
- The flutes are twisted, not ground, for superior performance, flexibility and strength.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

**Stainless Steel Pathfinder Instruments** feature a minimal taper and are used for negotiating difficult or calcified canals. Designed as a unique alternative to using #6 and #8 files, Pathfinder delivers more efficient cutting with less buckling, allowing you to use fewer files to negotiate difficult anatomy.

**Pkg of 6**

21mm	17023
25mm	17017

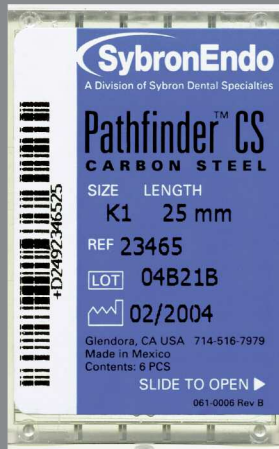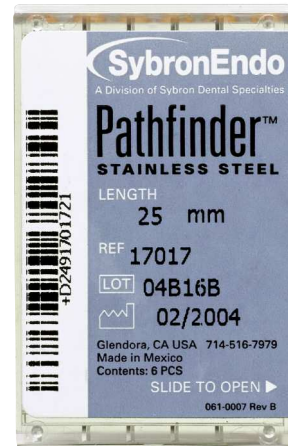

**Pathfinder CS Instruments** feature the same minimal taper as Stainless Steel Pathfinders but are made with carbon steel to help negotiate difficult or calcified canals faster and more efficiently than traditional stainless steel hand files.


**Pkg of 6**

Size	Handle	19mm	21mm	25mm
K1	Dark Brown	22935	22936	23465
K2	Orange	22937	22938	23467

**PATHFINDER CS**



- Made from carbon steel to create a stronger file.
- Perfect for breaking through calcification.
- Fewer files are needed – less buckling and more efficient cutting.
- K1 – between a #6 and #8 hand file.
- K2 – between a #8 and #10 hand file.
- The flutes are twisted, not ground, for superior performance, flexibility and strength.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

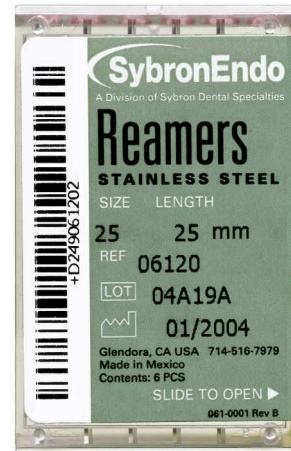

**SybronEndo Reamers** are high-quality, triangular-bladed, stainless steel instruments designed for flexibility with integrity.


Size	21mm	25mm
<b>Pkg of 6</b>		
08	14133	14134
10	06029	06110
15	06030	06111
20	06031	06112
25	06032	06120
30	06033	06126
35	06034	06127
40	06038	06131
45	06040	06132

**REAMERS**



- High-quality, triangular-bladed stainless steel files designed for better cutting efficiency.
- The flutes are twisted, not ground, for superior performance, flexibility and strength.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

**Frequent Filer Program**

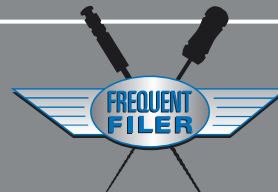


**As a reward for your loyalty,** you will accrue points for all of your file purchases. To learn more details about the Frequent Filer Program, contact your SybronEndo or dealer representative and **sign up today!**

# Strong, Safe and Versatile Stainless Steel Files



**OPEN THE PATH TO ENDODONTIC EXCELLENCE WITH K-FLEX**

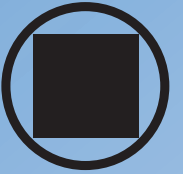


# SYBRONENDO FILES – THE CLEAR PATH TO PREPARATION SUCCESS



**K-FLEX FILES**

SybronEndo files and reamers are designed to cut quickly, efficiently and safely with maximized debris removal, even under difficult conditions. Highlighted by our innovative K-Flex Files, our comprehensive endo preparation line delivers a wide range of flexible yet rugged tools to easily handle your canal instrumentation needs.

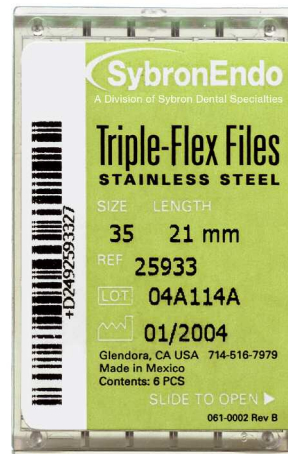

Size	21mm	25mm	30mm
<b>Pkg of 6</b>			
08	25927	25943	25959
10	25928	25944	25960
15	25929	25945	25961
20	25930	25946	25962
25	25931	25947	25963
30	25932	25948	25964
35	25933	25949	25965
40	25935	25951	25967
45	25936	25952	25968
50	25937	25953	25969
55	25938	25954	25970
60	25939	25955	25971
70	25940	25956	25972
80	25941	25957	25973
15-40 Assorted	821-3021	821-3025	821-3030
45-80 Assorted	821-3121	821-3125	821-3130

## TRIPLE-FLEX



- Triangular cross section for better cutting efficiency.
- Noncutting tip for safer canal instrumentation.
- Flutes are created by twisting the instrument rather than grinding, creating a more fracture-resistant shape and edge.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

**Triple-Flex Files** are triangular-bladed, stainless steel instruments designed for flexibility with integrity. The twisted flutes deliver the aggressive cutting and debris removal you demand without sacrificing the instrument's inherent stability. Unlike round blades that have flutes ground in, Triple-Flex Files are highly breakage-resistant, even in the most stressful situations. For performance and value you can count on, Triple-Flex Files are the preferred choice.

## K-FILES



- High-quality stainless steel square cross section for superior fracture resistance.
- Noncutting tip.
- The flutes are twisted, not ground, for superior performance, flexibility and strength.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

**K-Files** are designed to give the operator a smooth tactile sense inside the canal during instrumentation.




**K-Flex Files** (shown above) feature a unique rhomboid design, which provides a short axis for flexibility while maintaining cross-sectional strength. K-Flex is more flexible than most conventional files of comparable size, particularly in sizes #30 and above, which facilitates the negotiation and preparation of curved canals and minimizes potential for ledging and perforating.

Size	21mm	25mm	30mm
<b>Pkg of 6</b>			
08	15290	15122	15206
10	15294	15126	15210
15	15298	15130	15214
20	15302	15134	15218
25	15306	15138	15222
30	15310	15142	15226
35	15314	15146	15230
40	15318	15150	15234
45	15322	15154	15238
50	15326	15158	15242
55	15330	15162	15246
60	15334	15166	15250
70		15170	15254
80		15174	15258
15-40 Assorted	821-4021	821-4025	821-4030
45-80 Assorted	15368	821-4125	821-4130

## K-FLEX



- Diamond-shaped cross section yields improved sharpness, flexibility and cleaning.
- Alternating high and low flutes – low flutes improve debris removal and high flutes improve cutting ability.
- Instrument fluted to the tip provides better canal guidance and reduced ledging.
- Consistent instrument sizing and performance result in reduced instrument waste and expense.
- The flutes are twisted, not ground, for superior performance, flexibility and strength.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

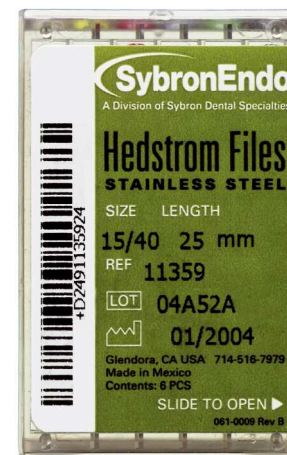


Size	21mm	25mm	30mm
<b>Pkg of 6</b>			
15	12766	11344	12814
20	12769	11345	12817
25	12772	11346	12820
30	12775	11347	12823
35	12778	11348	12826
40	12781	11349	12829
45	12784	11350	12832
50	12787	11351	12835
55	12790	11352	12838
60	12793	11353	12841
70	12796	11354	13278
15-40 Assorted	12809	11359	13291
45-80 Assorted	12810	11360	13292

## HEDSTROM



- High-quality stainless steel.
- Non-cutting tip.
- Stretch B Handle is designed with a serrated pattern that gives the clinician a comfortable, nonslip grip.

**Hedstrom Files** are ground from stainless steel wire. The flute depth increases with the taper, creating a sharper, more flexible instrument.

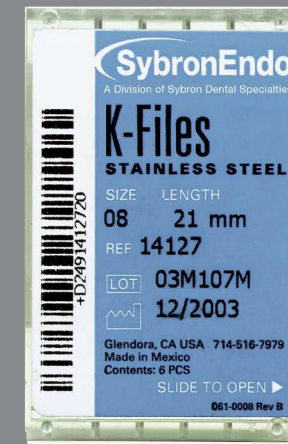

Size	21mm	25mm	30mm
<b>Pkg of 6</b>			
08	14127	14128	14129
10	06003	06062	06151
15	06004	06063	06152
20	06005	06064	06153
25	06014	06065	06155
30	06015	06066	06156
35	06016	06067	06157
40	06019	06069	06159
45	06020	06074	06160
50	06021	06075	06161
55	06022	06076	06162
60	06023	06077	06163
70	06024	06085	06166
80	06026	06099	06170
15-40 Assorted	821-8021	821-8025	821-8030
45-80 Assorted	821-8121	821-8125	821-8130